LUXCEAR Fornez 品質保証書

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LUXCEAR Fornez Warranty Certification

Please note that there will be no reissuing of this Warranty Certification, so keep it in a safe place so you may use it if necessary. If the product was purchased through a mail order or any other similar way, the date of delivery listed in the invoice will be considered the date of purchase, and the consignor the point of sale, so please keep the invoice in a safe place together with the Warranty.

* All stipulations of the Warranty shall be invalid if any of the spaces marked with an asterisk are not filled in, or if we cannot confirm the purchase with the retail store where you purchased your product.

Product Name	LUXCEAR Fornez		
Warranty period	1 year from the date of purchase		
*Customer Address	<input type="text"/>		
Customer Name	<input type="text"/>	Phone Number	(<input type="text"/>)
*Retailer Address	<input type="text"/>		
Store Name	<input type="text"/>	Phone Number	(<input type="text"/>)


★ Please understand that regarding the handling of customers' personal information, we may need to use the names, the addresses, and other personal information we are given by the customers for such purposes as consultations related to our company's products, repairs, and confirmations required therefor.
 ★ The Warranty refers to repairs offered free of charge within the period and conditions stated here. For that reason, the Warranty does not limit your legal rights. If you have any questions, please contact us directly at the inquiry desk.

Distributor: CONY Co., Ltd. 1-8-15, Azuchimachi, Chuo-ku, Osaka-shi, Osaka, 541-0052
 Inquiry Desk: 06-6262-0555 Open hours: 10:00 - 17:30 (excluding weekends and holidays)

Record of repairs

LUXCEAR Fornez

User Manual

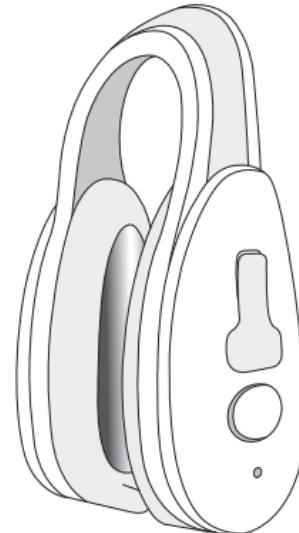
With warranty card

Thank you for purchasing this product.

- Please read this instruction manual carefully and only use the device in the correct manner.
- This document also serves as a warranty card. Ensure that it is stored safely and will not be lost.


1 Component Names
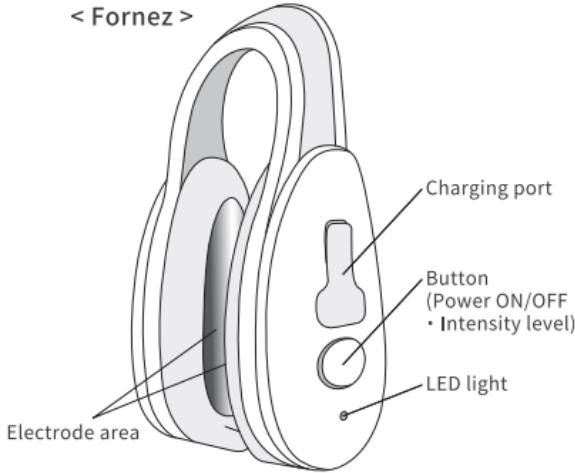
< Case >

< Charging cable >

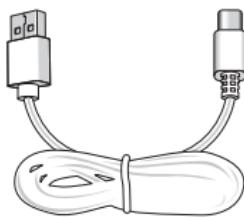2 Specifications

Fornez	Dimensions	W33 x H65 x H21 mm
	Weight	approx 19g
	Built-in battery	Lithium-ion rechargeable battery
	Charging time	approx 2 hours
	Continuous operating time	approx. 3 hours (when fully charged)
Charging cable	Length	approx 30cm

3 Safety precautions

- Please read the following safety precautions carefully before use and only use the device in the correct and proper manner.
- The precautions shown here are intended to ensure that you use the product safely and correctly, and to prevent harm or damage to the user or others.
- In order to clearly indicate the magnitude and degree of the potential harm or damage, the precautions are divided into "warning" and "caution" that are expected to occur if mishandled. All of these are important safety matters, so please follow them as directed.

Precautions for Users

Warning

Actions that may result in death (※ 1) or serious damages to health

Caution

The degree of harm or damage that is expected to result in serious injury to the user or to cause physical damage (※ 2) in the event of incorrect handling.

*1 Serious injuries are those that will have lasting after effects such as blindness, injury, burns (high or low temperature), electric shock, fractures, poisoning, etc. and those that require hospitalization or long-term hospital visits for treatment.

*2 Property damage refers to extended damage related to houses, household property, livestock, and pets.

Warning

- Never use with the following medical electrical devices:
Doing so may cause accidents or malfunctions of the medical electrical equipment, resulting in significant physical injury.
- **Implantable medical electrical devices such as pacemakers**
 - **Life-sustaining medical electrical equipment such as heart-lung machines**
 - **Wearable medical electrical devices such as electrocardiographs**

Warning

- Persons who have or have done the following should not use the device:
Doing so may result in accidents or skin or physical problems.
- Those who have or are suspected of having a heart disease (such as those who feel discomfort such as pain or palpitations in the heart)
 - Pregnant women and those who have just given birth
 - Those who are prohibited from exercising by a doctor
 - Those who have or are suspected of having venous thrombosis
 - Those who have consumed or are consuming alcohol
 - Those who have a fever (38°C or higher)

Precautions for Use

Caution

- Do not use in the following ways:
Doing so may result in an accident or injury.
- Do not allow use by small children or those who cannot understand the instruction manual by themselves.
Not doing so may result in accidents or skin or physical problems.
 - Keep this item out of reach of children.
Not doing so may result in injury or burns.
 - Never touch the electrode with your hands during use.
Strong stimulation may cause accidents or physical problems.
 - Do not disassemble, modify or repair the device on your own.
Doing so may result in an accident or injury.
 - Do not subject device to strong impacts, such as dropping or hitting against another object. Doing so may result in damage.
If you continue to use it as it is, there is a risk of accident or injury.
 - Do not allow any metals to come into contact with the electrodes.
Doing so may result in electric shock.
 - Do not touch the cable with wet hands.
Doing so may result in electric shock or burns.
 - Do not use with other electrical products. You may feel sick or feel unwell.

Caution

Please observe the following precautions before use.
Not following directions may result in accidents, injury,
electric shock, damage, or stains.

- Do not heat the device in a microwave oven and ensure it is placed away from fire.
- Do not use for anything other than its intended purpose or for anything other than what is listed in the instruction manual.
- Do not use on small children or pets.

Please do not allow small children and pets to play with this product.

- Do not allow others to touch this product while it is in use.

In particular, please do not allow small children and pets to touch it.

- Do not use in humid areas such as bathrooms.
- Do not apply excessive pressure to this product.

Caution


Please observe the following precautions before use.

- Should you experience pain or discomfort during or after use, or should you observe any abnormalities such as rough skin, rash, redness, itchiness, burning, or inflammation, discontinue use immediately and consult a doctor. If you continue to use it as is, it may worsen the symptoms.
- If this product emits a strange odor or noise, or if a foreign object such as water gets inside, turn off the power immediately. Not doing so may cause a short circuit, malfunction, or fire.
- Should you experience any abnormalities on your body or skin, discontinue use immediately and consult a doctor. Not doing so may cause an accident or poor physical condition.
- If you use it at an unnecessarily high intensity level or use it for a long time, muscle fatigue and/or a tingling sensation may be felt in the skin. Doing so may cause redness.
- When using this product, remove any accessories such as metal ones. If the accessory touches the electrode part, there is a risk of accident, injury, or damage.

4 How to use


Charge (Please charge the device after purchase and before use)

1. Open the charging cord port and connect the included USB C charging cable to the device.

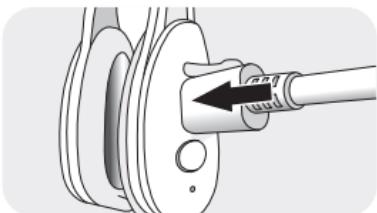
2. Connect the USB A end to an AC adapter*.

* This product does not come with an AC adapter. AC adapter must have a USB type A port (DC5V/500mA or higher output).
* Please charge in an environment with an outside temperature of 5°C~30°C.

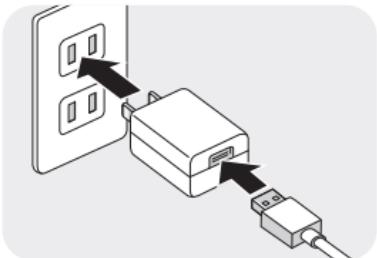3. Charging has begun when the red LED lights up. Charging is complete when the red LED goes out.
(About 2 hours to fully charge)

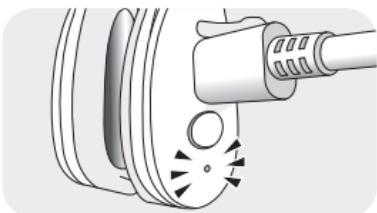Usage


1. Clean the nose thoroughly and remove oil, dirt, foundation,etc. from the are of use.

2. Press and hold the power button to turn it on. Press the button again to select a strength level. Please start from level 1. Please allow your body to adjust to the device before increasing the power intensity.


Level 1: Blue LED on

Level 2: Flashing blue LED

Level 3: Green LED on

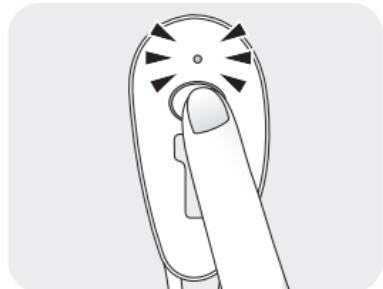
3. As shown in the figure, attach the wide side of the clip toward your face.
(Make sure the protrusion does not touch your nose.)

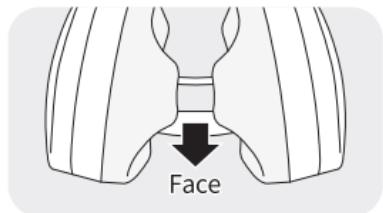
Usage

4. Attach it to a position that is secure and does easily slip.

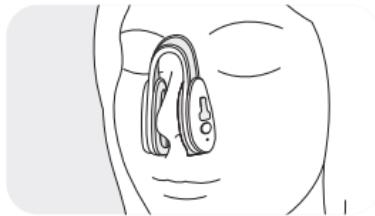5. It may be difficult to feel any stimulation at first, so please wait about 30 seconds after turning on the device.
(Different people feel the stimulation differently, but if you don't feel anything after 2~3 minutes, please increase the level.)

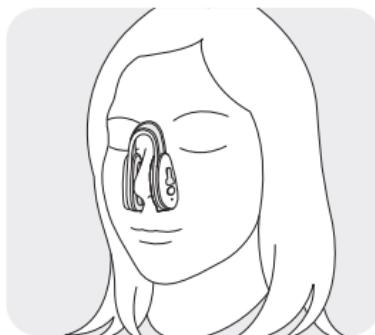6. Use at a comfortable strength.
Please use for 10~15 minutes at a time.
The power turns off automatically after 15 minutes.


How to Remove the Batteries

* Do not remove the battery for any purpose other than for disposing this product.

1. Make sure the power is off.
2. Insert a flat-blade screwdriver and remove the cover using some leverage.
4. Remove the battery by pulling it out of the connector.

3. Remove the screw that is fastened to the board.
5. Fasten the end of the connector with insulating tape (cellophane tape, etc.).

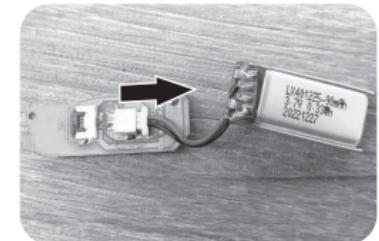Troubleshooting (Q&A)

Q There is no stimulation at all; Is it a malfunction?

A: Please thoroughly remove dirt, sebum, and foundation from your nose before use.
It is also recommended to use skin toner to improve the conduction of electricity.

Q Do you have to remove the silicone part of the device to use it?

A: No, please use it as it is without removing any parts.

Q It's hard to understand which way I should wear the device.

A: The charging port should be on the left when you put on the device.

Q I feel some pain in my skin, should I continue to use it?

A: Please refrain from using it until the pain disappears.
Then, after your skin condition has settled, use it only for a short time at the lowest level.

5 Device Care

- ① After use, gently wipe off the electrode with a dry towel or tissue.
- ② Please store it in the provided case.

6 How to dispose of this product

This product has a built-in lithium-ion rechargeable battery.
When disposing of it, please follow the ordinances of your local government.
For details of the ordinance, please contact your local government.

7 Warranty and after-sales service

■ Warranty Card

Please read the warranty card carefully and keep it in a safe place together with the proof of purchase (receipt, etc.).
The warranty period is one year from the date of purchase.

■ Inquiries about this product

If you have any questions about this product, please contact us at:

Distributor: CONY Co., Ltd. 1-8-15, Azuchimachi, Chuo-ku, Osaka-shi, Osaka, 541-0052

Inquiry Desk: 06-6262-0555 Open hours: 10:00 - 17:30
(excluding weekends and holidays) (Service in Japanese only)